

## **Power Technical Training**

Including Specialty Controls and IR

Ken Kalinowski, Senior Service Technician ken.kalinowski@sunmed.com















## Agenda & Objectives

### Course Outline

- Course Outline

  Intro Tech Webste Overview

  Mid Wheel Drive & Rear Wheel Drive Chair Quick Overview
  Electronics

  VR2 Standard Electronics

  R-Net 90 Amp PM-120 with All options

  Programming
  R-Net including OBP On Board Programming
  PC Programming
  Specially Controls
  Setup and Operation
  Blue Took Stdup and Operation
  Input / Output Module
   Usage & Aquatiment
  IR Setup and Operation
   Delete and Remane
  Assignable Buttons
   Setting Second Fraction
  In Public Second Fraction
   Setting Assignation and Prover

  Seating ASAPII Seating
  Seating Ontrols

  Seating ASAPII Seating
  Seating Controls



## The People Behind the Power







STEPS





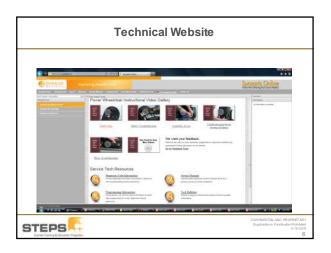

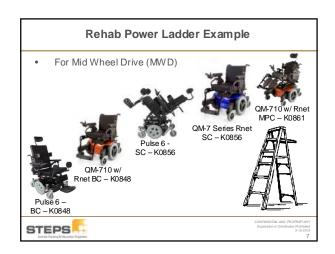





## Power: Simple and Reliable

Is the manufacturer focused?

Do they specialize in Group 3 Power with varying models?

- Proven Reliability in Components
- Motors
- Electronics
- Tilt Module
- Made in the USA



## **Electronic Selection**

- Standard Level Package
  - Penny & Giles VR-2
    - Programmable
      - With hand held or PC programmer
      - 1 profile with 5 speed increments
      - 5 profiles with 1 programmed speed
    - 90 amp controller
      - Uses the same parameters as Qtronix
      - Provides activation of 2 power options through the joystick
    - Controller to be used with a joystick only
      - Not expandable.





## Standard Electronics Package

- VR2 Electronics Operation
- 4 or 6 Button Joystick:
  - This joystick has a GREEN On-Olf switch with RED horn Top display gives indication of both battery capacity and Diagnostic Codes. The center indicators are 5 individual YELLOW LED lights that indicate

  - The determination of which does in the control of the state indicate


     Either Speed (if module is set to "0" profile)
     Or inshidutal lights indicating separate drive profiles up to 5.

    The operation of those indicators is accomplished by pressing the right Speed Drive select to increase or the left Speed Drive select to decrease the speed or the left Speed Drive select to decrease the speed of the left Speed Drive select to decrease the speed of the left Speed Drive select to decrease the speed of the Actuator or Seating buttons indicate via a RED LED that it is active and a Forward and Reverse action to operate the actuator.

    The 4-trin connector connects to a compranication 4-trin
  - The 4-pin connector connects to a communication 4-pin connector that connects back to the Motor Control.





STEPS

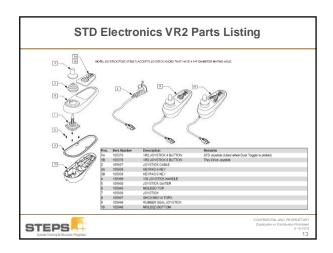







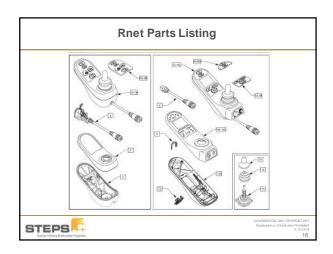


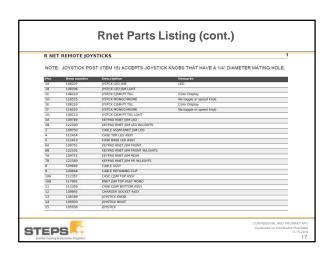

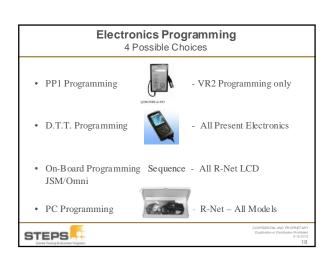


## The DTT (Diagnostics Test Tool) has USB port. Files can easily be transferred using a USB drive. - Holds up to 1000 profiles The Diagnostic Test Tool has been designed specifically to assist technicians, enginees and specialist. The DTI drives uses to access, program and show Po Dieve fectoricially specific files in a conversired and effective way. The Diagnostic Test Tool will drive uses to. - Program Riu Dieves Technology controllers. - Read and save program files. - Read and save diagnostic logs. - Read and save diagnostic logs. - Read and save diagnostic logs. - Save and delete program files boothy. - Save and delete program files boothy. - Save and delete program files boothy. - Save and delete program files toothy. - Save and delete program f

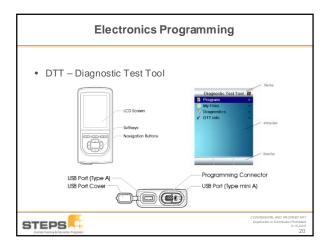

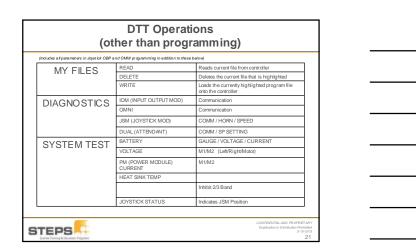


# DTT Operations Continued (other than programming) (Includes all parameters in Jugsteck OBP and OLMS programming in addition to these below) SYSTEM LOGS | IOM / OMNI / JSM / DUAL / PM | Lists faults and codes of connected modules | MY LOGS | Save the System Log file into My Logs | SYSTEM | How many hours has it been driven | TIMERS | DTT INFO | Contains information pertaining to the DTT, including the software version.



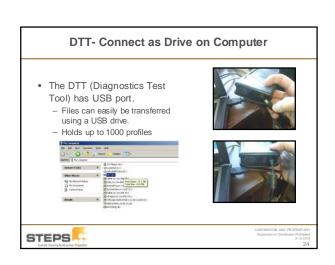



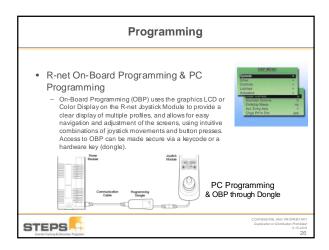



## On Board Programming Mode Key Sequence


- If the keycode method is set, then the following button sequence will allow entry to OBP:
  - Note, a dongle will also allow OBP access if this method is set.
  - Hold down the Horn button and then hold down the On/Off button until there is a short bleep.
  - 2. Power-up sound will occur prior to this bleep.
  - 3. Release the Horn button, but continue to hold down the On/Off button until there is a further short bleep.
  - 4. Release the On/Off button, there will now be a longer bleep and OBP mode will be entered.



CONFIDENTIAL AND PROPRIET AS Duplication or Distribution Prohibits 3-15-20

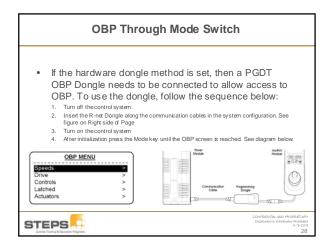

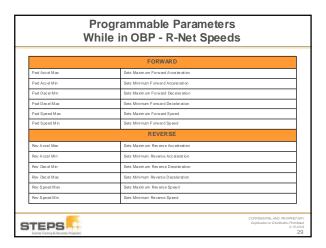

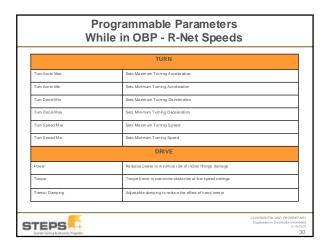


	GLOBAL CONTROLS	
Steer Correct	Adjusts the P M output's to compensate for mis-matched motors	
Sound V olume	Sets volume of audible feedback from JSM.	
Endst op Beep	Sets whether there is a bleep when a seat axis reaches endstop	
Act. Entry Axis	Sets the default axis when seat control mode is entered	
Chge Prf in Dr v	Sets whether profile changes are permissible while driving	
Speed Adjust	Sets whether the speed buttons on the JSM are active	
Speed A dj in Drv	Sets whether speed setting changes are permissible while driving	
Moment ary Screens	Sets Whether large screens appear at profile and speed changes	
Rev Driving Alarm	Sets if the reverse driving alarm is active	
Lock Function Enabled	Sets how the lock function is activated	
Display Speed	Sets how the speedometer is displayed (miles or kilometers per hour)	
Max Display Speed	Sets the operation of the graphical speed display	
Power – Up Mode	Sets the Mode that will be active when the system is powered up	
Display Speed	Displays the Speed – MPH – KM/H – or Off are choices	
Max Display Speed	Displays Max Range of Speed - Sets 0.1 to 20 - set at motor speed	


Controls				
Controls				
	PROFILED CONTROLS			
Sleep Timer	Sets the time of inactivity before the system goes to sleep			
Change Mode in Drive	Sets whether mode changes are permissible while driving			
Back ground	Sets the default back ground for each profile joystick			
	JOYSTICK			
Active Throw Bds psystick throw via joystick movements Throw Detail Sets joystick throw via programming				
				Active Orientation
Orientation Detail	station Detail Sets joystick orient ation via programming			
Deadband	Sets the joystick deadband (size of neutral position)			
		CONFIDENTIAL AND PROPRIET  Duplication or Distribution Profi		


STANDBY				
Switch to Standby Sets whether an external button can be used to enter Standby Mode				
Standby Time	Sets the time of inactivity before Standby Mode is entered			
Mod e Se le ct	Sets whether of her modes can be selected from Standby Mode			
Remote Select	Sets whether a profile can be selected from Standby Mode	┪		
	LATCHED			
Drive	Selects latched drive operation			
Actuators	Selects latched a ctuator operation	Selects latched a ctuator operation		
Timeout	Sets the time out period for latched operation	٦		
Timeout Bleep	Sets whether a bleep occurs as the timeout period approaches	٦		


	SM ACTUATORS				
Away Speed	Sets actuator axis speed in the up direction				
Hom e Speed	Sets actuator axis speed in the down direction				
	ISM ACTU ATORS				
Away Speed	Sets actuator speed in the up direction				
Hom e Speed	Sets actuator speed in the down direction				
Acceler ation	Sets actuator acceleration				
Dec el erat ion	Sets actuator deceleration				

Controls					
	PROFILES				
Profile E nable	Sets which profiles are available to the user				
MODES					
Mod es	Modes Sets which Modes are available in each profile				
INPUT DEVICES					
ID Type	Sets which Input Device controls a Profile				
SYSTEM					
Diagnostics	Accesses the system's error log				
Timers Accesses the system's run timer					
Joystick Calibration	tallity to copy one profile to another  Next one factory preset values "  Rest one factory preset values "				
Copy Profiles					
Load Presets					
Adjust Contrast					
TEPS		CONFIDENTIAL AND PROPRIE Duplication or Distribution Pro 3-1			

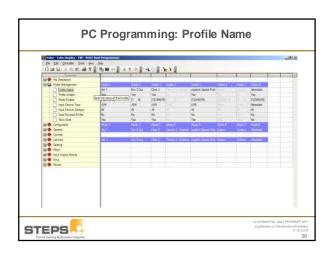

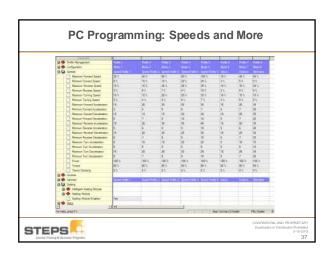

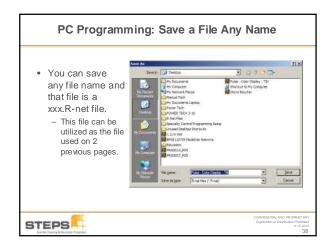


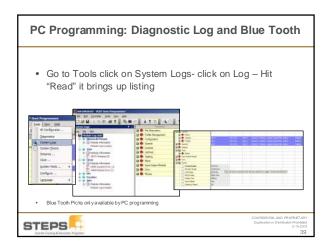



## Specialty Controls COMPERITM ARE PROPERLY AND PROPERLY PROPERLY AND PROPERLY PROPERLY AND PROPERLY PR

## IR (Infra Red) control is standard with all OMNI control packages The wheelchair user can use any Input Device connected to an OMNI control for operating IR controlled devices such as TVs and DVDs. Up to 100 programmable commands Bluetooth communications for PC Mouse Control with any Bluetooth compatible computer platform. Bluetooth control is the standard platform used for the wireless mouse. Operates from up to 10 meters away − no need for line of sight. CONCENTRA NO PROPRET ANY CONCENTRA NO

# Specialty Control: Section Features General R-net interface for specialty controls Large, hi-resolution color display Dual channel inputs User selectable display - indoors/outdoors Customisable user menu Large, easy-to-read icons Programmable user text On-board configuration/programming Use standalone or with other R-net Input Devices Can control all powerchair functions

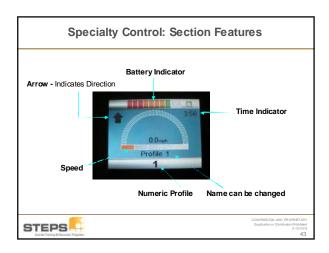

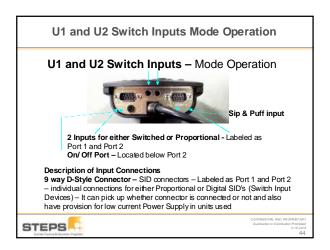


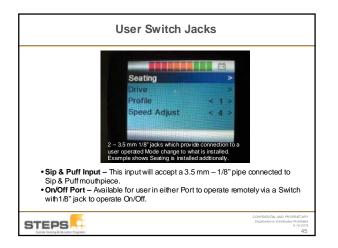







## LCD Screen A full color, backlit LCD screen that can show configuration details and operating information to both dealer and user. ON/ Off Button This button provides a complete Power Down of control system electronics. In addition to the button mounted there is a facility on the inputs for an optional, externally mounted switch that can be operated by the user for on/off usage. Mode Button The mode button allows for the user to change between available programmed Modes. Profile Button The profile button allows for user to change between available programmed Profiles or Drives. Navigation Buttons The array of 4 buttons allow Dealer to navigate between OBP screens.

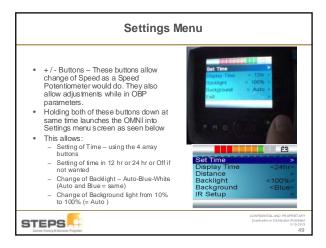


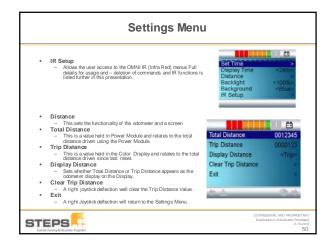


## 

**Program OMNI To Match SID** 

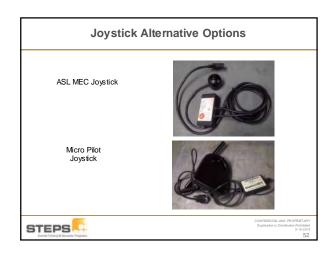

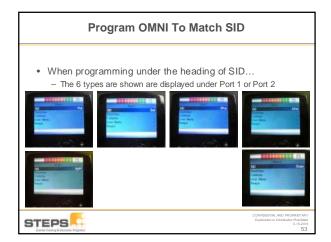

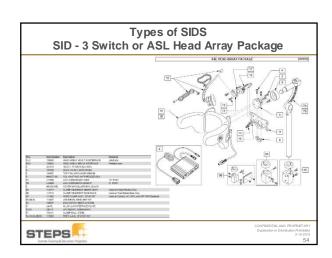

## R-Net Specialty Control Set-up Guide

- 3 Switch Head Array
  - (Please use these instructions as a guide only. Due to multiple options, it is important to select the correct set-up to meet the individual needs of the user.)

     Enter On-Board Programming (OBP) Thru the OMNI Display
  - - Connect up Dealer dongle to cable assembly
       Reset the Power and Press the Mode Button.

    - . Press the Mode button until OBP comes up on Screen follow details below
- Step 1 Setting up the Motor Controller Drive Profile(s) for the appropriate Specialty Input Device (SID)

   Scroll down using the directional buttons on the OMNI display to Input Devices. Use the Right navigation button to select.

  - The input device either needs to be "Uni" for universal input, or "Omn" for OMNI Display. Use the + & buttons to select the desired SID.

    Left Navigation button back to Main OBP menu. This will save your changes.



## R-Net Specialty Control Set-up Guide

- Step 2 Set Up From OMNI
   From the main OBP menu, highlight 'OMNI' and select by pressing the right directional button on the OMNI display.
  - right directional button on the OwiNi display.

    Scroll down to "Profiles" using the directional buttons, and select by pressing the right navigation button.

    There are two Specialty Control Ports on this OMNI display, when setting up a specialty Control we need to identify which port will be used.

    Set all profiles you would like to control with the Specialty control device to Port1. Left navigation button twice to save your changes.

  - Use the down navigation button to highlight "Port 1", and select using the right navigation button.

    "SID" is the first programming option in this menu. Use the + & signs on the OMNI to select "3Swi".
  - Scroll down to Switches and select using right navigation button.
     Confirm "User Switch" is set to "NO" (normally Open). (If it is not, change to "NO" using the OMNI's "+ & -" buttons)

  - Left navigation button twice Back to the OMNI OBP menu.



## R-Net Specialty Control Set-Up Guide

- · Reverse Options:
  - Double Click Time programmable from 0-2.5 Seconds in 1.0 increments

     A single operation of the user switch toggles directions.
  - A double operation of the user switch will enter into user menu
     Forward / Reverse Auto Toggle Programmable On or Off

  - Direction change will occur if a forward/reverse command is operated and released (rear pad of headrest).
  - Further operation of the of the rear head pad will result in drive in the newly selected direction.
  - Switch Medium Programmable from 0-5 in .25 increments

     A directional change occurs at the instant the user switch is operated.

    - If the user switch is released and a forward/reverse is given, then drive will commence in the newly selected direction.

      Double Click Time must be set to 0

      - Make sure your Switch medium & switch long settings are not the same




## R-Net Specialty Control Set-Up Guide

- Automatic Menu Scan:
  - If you want to program an Automatic Menu Scan so when the client selects mode for a drive profile, the User Menu with automatically scan thru the options..
  - this can be done in Std. programming. Go to OMNI / Port 1 / Controls / Menu Scan Rate.
    - Menu Scan Rate is programmable in % second intervals from 0 to 10 seconds. (This function is Off when set to 0)
      - If your mode button is built into your switched input device you may also need to turn off the User switch detect.
    - From the Main Menu in OBP... Highlight "OMNI" and select using the right navigation button.
      - Highlight "Port 1" and select.
      - Highlight "Switches" and select.
      - Highlight "Switch Detect" and using the + / buttons change to "OFF".





STEPS

## **Programming Functions for OMNI Control Global Settings** Puff Threshold Put 1 Troeshold The sace is the rehelded point between soft and hard puts – 10 – 10.0 – 1.0 increments The sace is the rehelded point between soft and hard puts – 10.0 – 1.0 increments The sace is the treathed point between soft and hard puts – 10.0 – 1.0 increments Double Bard The sace is the size of the spipul material are a how much pressure the user must apply with either a sign or a putf before the brakes Double Click The presentance or the spipul material are a how much pressure must apply with either a sign or a putf before the brakes To bouble Click The presentance operations made with the presentance or the presentance of the spipul material are a deather held or command to make the spipul material are spipul materials. Sets the period of the which help presumate commands must be described. Sets the period of the which help presumate commands must be described.

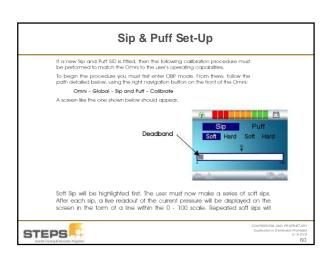


# Sip & Puff Set-Up • The user must now make a series of hard sips to produce a new Hard Sip band. Ideally, there should be as much difference as possible between the Soft Sip and Hard Sip values. • The default for Threshold is 50 – once values are saved it is mentioned for optimizing to set Threshold not the mid point between soft and hard. This is done with the +3 and – buttons on the front of the OMNI. • Now repeat the sequence of instructions as above for both soft and hard puffs including the Threshold Setting. Once calibration button on the front of the OMNII to return to the Calibrate option in the Sip and Puff branch of OBP.

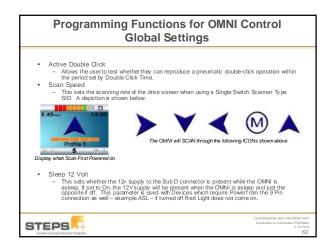


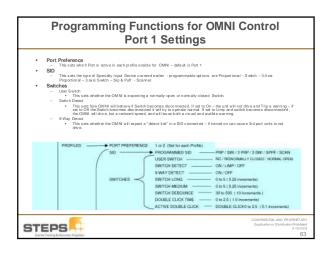

## Programming Functions for OMNI Control Port 1 Settings Switches - Switch Long \* This sets the period of time the User Switch must be operated in order to put the OMNI into the sleep mode. 0.5 sec to 5.0 sec in steps of 0.25s. Switch Medium \* The parameter is intended for 3 Axis Proportional and Switch Type SIDS only - direction change occurs at the instant the User Switch is operated, provided the parameter Double-Citck has been set to 0.1 User switch is referse and a Forward Reviews command is entered then drive will continue in the newly selected direction. If the User Switch is not released and is held for a time period set by the operation. 0.5 sec to 5.0 sec in steps of 0.25s. Switch Debounce \* This sets the debounce time for the Userswitch Debounce is period of time the switch is operated continuously before a new condition can be registered. 30 – 500 – increment 10 Double Click \* This sets the period of time during which two User Switch clicks must be detected. In order to interpreted as a double-Click. 0 – 2.5 – 0.1 increments Active Double Click \* This parameter allows user to set and test the Double Click Speed. Note: This parameter is only available in OBP.

STEPS

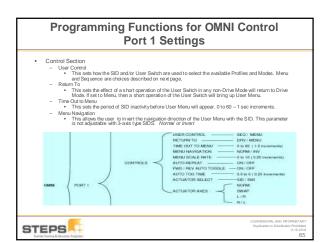

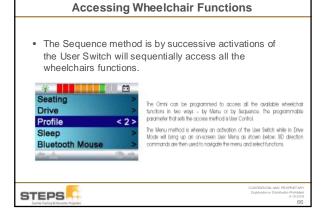

Controls - User Control

## **Programming Functions for OMNI Control Port 1 Settings** Control Section Menu Scan Rate This sets the scanning rate for the User Menu, Lights Menu, Menu and Settings Menu. O to 10 -- 0.25 increments - set to 0 no scanning and default. Auto-Repeat This sets whether there is an auto-repeat function on the direction commands for various SIDs. Auto-repeat means that successive commands will be applied if the switch is cortinuously operated, which is useful when navigating the User Menu. Not desirable for head control applications. On — Off Fwd/Rev Auto Toggle SIMEN ALLO Toggle This sets another method of changing direct with 3-Axis. SIDs. This has no effect with any other type of SID – On or Crt. If set to On, than an operation of the SID's FwdRev command can be used to change the selected drive direction. This is achieved by operating and releasing the FwdRev command within a 2 second period. This will change direction of drive, then the FwdRev commandmust be operated again within a further 2 seconds. If this 2 second period passes without a FwdRev command, the selected drive thereaften will revert to original beating. Note: This will not work with Latch and will result in a Bad Setting Error appearing. Actuator Selection This sets whether actuator is selected via SID commands or User Switch



## **Programming Functions for OMNI Control Port 1 Settings**

- Control Section
   Actuator /xis
   This sets which SID direction commands are used to select available axes. Not available for Single Switch Scanner. Settings for Actuator Axis are as follows:

  - Normal

    Left and Right commands will select available actuator axis, and forward
    and reverse commands will move actuator(s).

    Swap

    Forward and Reverse commands are used to select available axis, and left
    and right commands will move the actuator(s).

  - A left command will select the available actuator axis and right command will move the selected actuator(s). The reversal of operation (toggling) is change via a short operation by a right command.
  - - ght/Left

      A right command will select the available actuator axis and left command will move the selected actuator(s). The reversal of operation (toggling) is change via a short operation by a right command.





## Using the OMNI With Two SIDS

- There are 5 programming Steps to Follow:
  - 1. Decide which Profiles are associated with OMNI whenever deciding on SIDs.
  - 2. Decide which Port to select for each of the SIDS a selection of Port 1 or Port 2.
  - 3. Set the Parameter Profile Enable only required profiles.
  - 4. Ensure the Parameter Mode Enable is set to enable the required Modes for Both SIDS.
  - 5. Decide upon appropriate names for each of the Profiles and set the Profile Name accordingly using the PC Programming.



# OMNI Programming (Avail. When Accessing OBP through OMNI or DTT) • Includes all parameters in Joystick OBP programming in addition to these below: | Parameters | Parameter

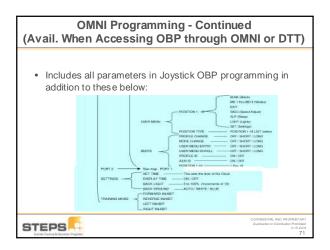

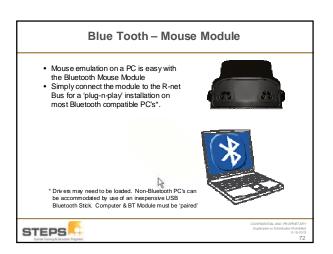

## Blue Tooth - Mouse Module **R-Net Only**

- You will need
  Power Module
  Bluetooth Module
  Input Device Module e.g. JSM or OMNI with a specialized input device like a Penta switch.
  A standard JSM will be assumed for the rest of this example.

  - When you connect the system for the first time it will go through configuration "handshaking" and ask to be turned off and on again. When this is completed erter Bluetooth mode by pressing the 'Mode' button on the JSM. Bluetooth mode is represented on the LCD by a laptop symbol displaying the Bluetooth logo.
  - Slowly deflect the joystick in a full forward direction and hold for 10 seconds, at which point a beep will be heard. (N.B. entering Bluetoothmode and deflecting the joystick must occur within 10 seconds of the system being furned not to work correctly).
  - 4. After the beep, deflect the joystick in full reverse until it beeps again (about 10 seconds).

  - 6. R-net Bluetooth is now in Discovery mode, so can be picked up as a Bluetooth device by the



## Blue Tooth - Mouse Module **PC Setup**


- On your PC, in Control Panel go to 'My Bluetooth Places'.
- · When you open 'My Bluetooth Places' click 'Add'.
- Tick the box 'My device is set up and ready to be found'.
- This will tell the PC to pick up Bluetooth devices in range.
- A device called 'PG Drives Bluetooth Mouse' should appear in the window.
- Click on the 'PG Drives Bluetooth Mouse'. This brings up a window called Bluetooth wizard with 4 options relating to how to connect using passkeys.
- · Select option 3, and enter the Pairing code '0000'.
- This pairs the R-net Bluetooth mouse with the PC.





## Blue Tooth - Mouse Module **Programming on PC**


- There is programmability for users that cannot operate the Speed Up / Down buttons. This is accessible via the Bluetooth branch of the PC programmer. (Shown on Next Page)
- The parameter is called 'Nudge'. A Nudge is defined as a quick deflection of the joystick or Specialty Control in one direction.
- There are 4 Nudge functions, each of which can be programmed to act as Scroll Up, Scroll Down, Left Click, and Right Click:
  - 1. Forward
  - 2. Reverse
  - 3. Left
  - 4. Right
- There is "New" Pointer Speed and Acceleration that allow changes to speed of Pointer and the Acceleration across the screen. There is "New" External Switches that are available for clients that cannot get the Nudge function but can hit a Switch.



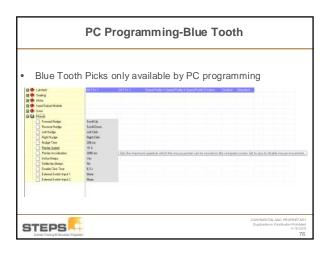

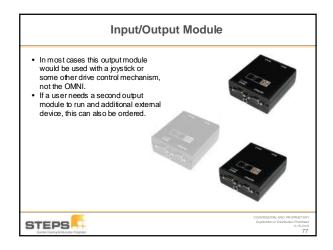

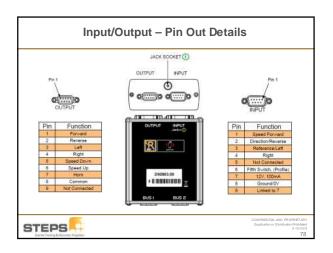

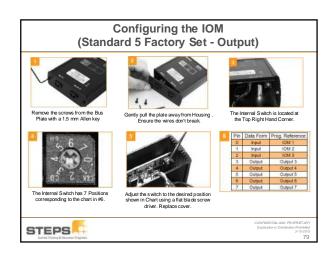


## **Input Setup Configuring for Input Device**


- To fit a Head Array Switch unit, without a Rebus Connector, into a control System:
  - Set the internal switch to Position 0.

  - 2. Connect the IOM to the R-Net control System.
    3. Connect the Head Array to the IOM 9 Pin Sub D Input Connector.
  - Turn the Control System on using the Head Array On/Off Control, or an On/Off switch connected to the IOM's On/Off jack socket.
- Details: The PM or Power Module is factory set to give control to the Input Device that has been used to Power it up. It is therefore possible to connect multiple input devices via IOM's. It is also possible to assign specific Input Devices to specific Profiles




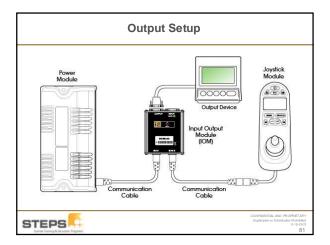



## Output Setup Configuring Output Device & Specs Output Functions: To fit an Environment Control Module, without Bus Bar, into a control system. Set the internal switch to Position 5 Connect the IOM to the R-Net Cortrol System. Turn the Cortrol System on and use the Mode Button to move to Mode 5. The PM (power Module) – is factory set to accept IOM Output Devices in Modes 5, 6, and 7 when the internal switch is set to 5, 6, 8.7 respectively. It is possible to connect up to 3 IOM as Output devices. You may have to adjust the programming parameters to ensure Injut device is selected accordingly. Diagonal Control Only one output can be switched at one time. The output direction is joystick xy data is shown on the side.

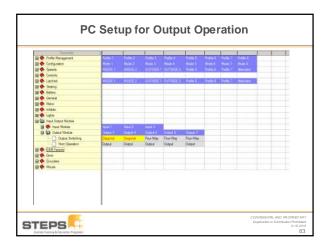




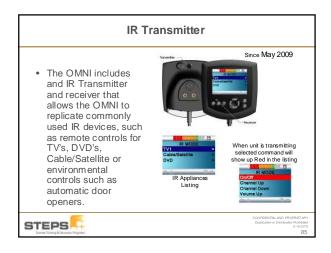

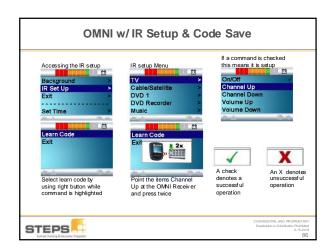

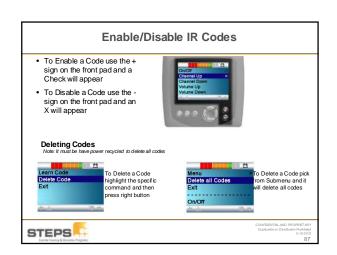

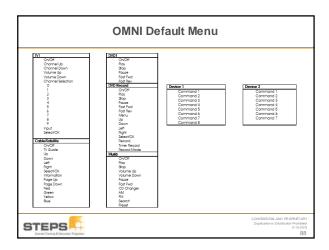

## **IR Configuration Tool**

- The IR configuration tool is a PC based application that allows the user to
  - Read and Write IR menus from and to an R-net Control System
  - Create IR menus.
  - Change IR menus.
  - Save IR menus.




CONFIDENTIAL AND PROPRIET AF Duplication or Distribution Prohibit 3-15-20


## **IR Configuration Tool**


- The IR configuration tool is a PC based application that allows the user to
  - $\,-\,$  Read and Write IR menus from and to an R-net Control System
  - Create IR menus.
  - Change IR menus.
  - Save IR menus.





CONFIDENTIAL AND PROPRIET A Duplication or Distribution Prohibi 3-15-21

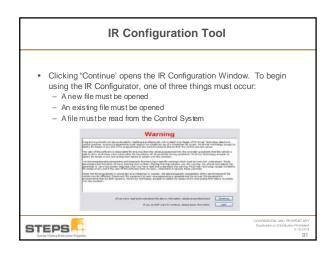


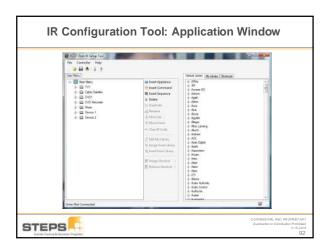

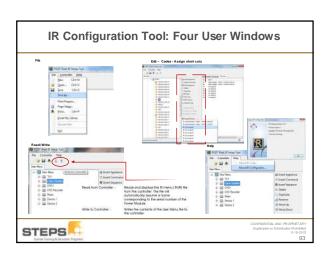

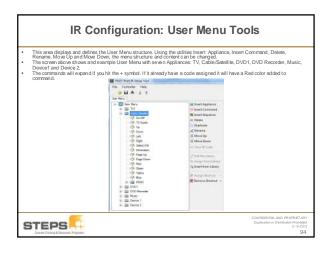

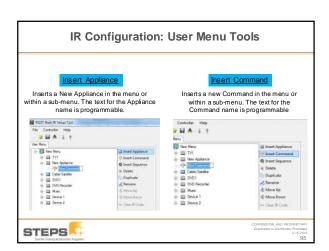

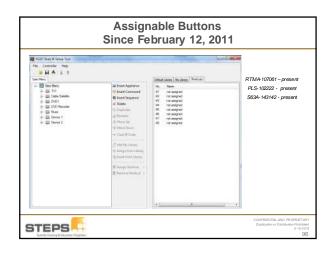

## **Assignable Buttons - Overview**

- The Assignable Buttons feature is a proprietary feature owned by Sunrise Medical.
   It enables the re-assignment of buttons and jacks on joysticks and specialty control modules.
- modules.
  When using Assignable Buttons, selected buttons and jacks can support up to two features.
  The second feature is controlled by time.
- If a button or jack is held closed for longer than the Second Function Time, the assigned second function will be activated.
- The color joystick chosen by Sunrise NA uses speed pots which disable the Speed Up and Speed Down buttons.
- Using the Assignable Buttons feature users will now be able to use the Speed Up and Speed Down buttons for a dealer settable action. By default he second function of the Speed buttons will be Tilt up and Tilt down
- Special tribes second tribution of the Special buttons will be filled and fill trespectively.

  Color Joystick now supports stereo-jacks in both positions. This allows the assignment of two functions to the Profile/Mode jack.
- The functions of the On/Off jack are not changed and are not assignable.



## Assignable Buttons: Breaking News

- A late breaking feature for Assignable buttons is that the second function time is dealer settable.
  - The Dealer PC programmer upgrade is required to correctly set the assignable buttons features.
    - This can be downloaded from a link to PG Drives from the Sunrise website from the product pages -> PG Drives Technology



## Assignable Buttons: Hardware Overview

- The following input modules support Assignable Buttons
  - LED JSM

  - Color JSM
     OMNI Specialty Control Display
     Bluetooth module
- · Color joysticks with stereo jack.
  - This is a new hardware feature being introduced to NA along with Assignable Buttons.
  - Stereo Jacks have been installed in color joysticks delivered to Suniise beginning in June 2010
- Monochrome joysticks do not support Assignable Buttons.
   Only the LED and Color joysticks have implemented the Assignable Buttons feature. Monochrome joysticks do not currently contain stereo jacks.
- Assignable Buttons features are not available with on-board programming (OBP) or with the Diagnostic Test Tool (DTT)



## Assignable Buttons: Software Overview

- Controller Personality (rnd) file.
  - This file is loaded into Sunrise controllers during the end of line process.
  - It is also necessary to be activated within the PC Programmer (either Dealer or OEM) to display or modify the Assignable Buttons features.
  - The Controller Personality file controls access levels of parameters and features.
- Assignable Buttons features will be enabled on all chairs.
  - Second function time will be set to 1 second. This is an OEM only
  - parameter.

    Dealer kits will continue to contain the standard version of the Dealer PC Programmer

    Sunrise Medical has a link to the Sunrise Dealer PC Programmer.


  - Latching actuators are not required at this time but may be implemented at a later date. The next page gives details on how to download the latest dongle software.



## Website Link to Upgrade: Dongle Software SUNRISE QUICKIE Live without limits the class Printy Compile Aurent Cost of Compile and Minimay Notes: Bent Teachers: Typing Regulations. STEPS






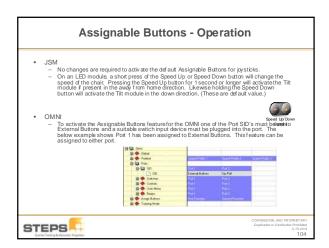




## The button functions can also be assigned a toggle and/or latched capability Rule of thumb!! Don't put anything that will run for any length of time (like and actuator) on the first function. Otherwise a may take 20 hits to tilt a seat back. Example of direct functions Example of toggle Short hit (less than 1 second) switches modes – direct function Long hit (more than 1 second) tilts seat back – direct function Long hit (more than 1 second) tilts seat back – direct function STEPS STEPS Convertible And or MORES FMY Children in All and or

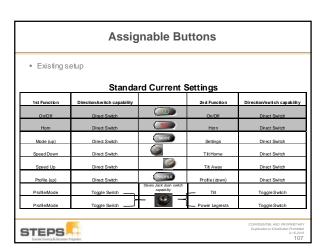



	Assignable Buttons						
	Beginning <u>Sept - 2011</u> these second function settings will be based on the chair configuration						
١.	Tilt Settings						
	1st Function	Direction/switch capability		2nd Function	Direction/switch capability		
ΙĮ	On/Off	Direct Switch		On/Off	Direct Switch		
[	Hom	Direct Switch		Hom	Direct Switch		
	Mode (up)	Direct Switch	MODE	Settings	Direct Switch		
	Speed Down	Direct Switch		Tilt Home	Direct Switch		
	Speed Up	Direct Switch		Tilt Away	Direct Switch		
	Profile (up)	Direct Switch	PROFILE	Profile (down)	Direct Switch		
	Profile/Mode	Toggle Switch —	Storeo Jack dual- switch capability	— Tilt Hame	Toggle Switch		
	Profile/Mode	Toggle Switch	_	Tilt Away	Toggle Switch		
UJ.	STEPS  CONFESTML NO PROPER FAVY  Deplication of Database Player  Local Confest Annual Principal  Confest Annual Principal						


Power	Legrest Se	ettings	
Direction/switch capability		2nd Function	Direct io n/switch capabilit
Direct Switch		On/Off	Direct Switch
Direct Switch		Horn	Direct Switch
Direct Switch	MODE	Settings	Direct Switch
Direct Switch		Left/Extend Leg	Toggle Switch
Direct Switch		Right/Extend Leg	Toggle Switch
Direct Switch	PROFILE	Both Legrests	Toggle Switch
Toggle Switch	Stereo Jack: dual- switch capability	- LeftExtend Leg	Toggle Switch
Toggle Switch		Right/Extend Leg	Toggle Switch
	Direction/switch capability  Direct Switch  Direct Switch	Direction/twitch capability  Direct Switch  Sizero Auth Calair switch	Direct Switch  Direct


Assignable Buttons							
1st Function	Tilt & Po	wer Legrest	Settings 2nd Function	Direction/switch capability			
On/Off	Direct Switch		On/Off	Direct Switch			
Horn	Direct Switch		Hom	Direct Switch			
Mode (up)	Direct Switch	MODE	Settings	Direct Switch			
Speed Down	Direct Switch	<b>(</b>	Left/Extend Leg	Toggle Switch			
Speed Up	Direct Switch		Right/Extend Leg	Toggle Switch			
Profile (up)	Direct Switch	PROFILE	Both Legrests	Toggle Switch			
Profile/Mode	Toggle Switch —	Stereo Jack dual - switch capability	- Left/Extend Leg	Toggle Switch			
Profile/Mode	Toggle Switch		Right/Extend Leg	Toggle Switch			
				•			
TEDS				CONFIDENTIAL AND PROPRI Duplication or Distribution Pro			

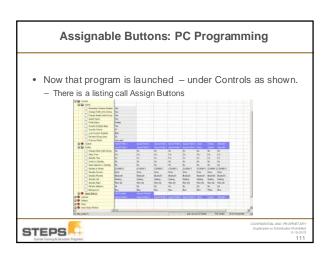


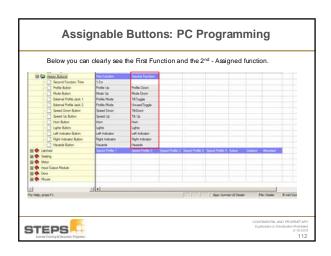

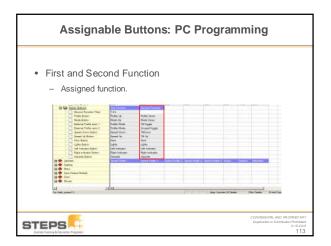

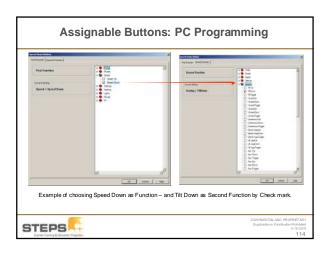


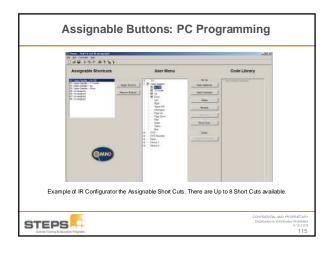


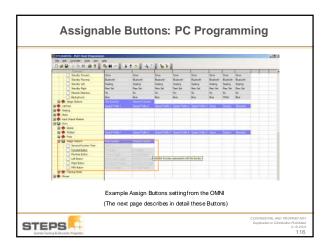

### To suit specialist user preferences, up to five of a Specialty Input Device's buttons can be programmed to have differing first and second functions. First Function Operates on a short depression and release of the button, while the second function operates if the button is held depressed for a period of time. Second Function (The Parameter) Time Sets the period of time the button must be depressed for its second function to operate. Note, if the first and second functions are set to be the same, then the function will operate on the depression of the button. The buttons that can be programmed with different functionality include: Forward Reverse Left Riight Fifth To set the button functions, fiere are programmable parameters with the same names as the buttons and each has a First Function and a Second Function setting.

# Enables Stereo Jack on Color JSM. Does nothing if applied to a unit without a Stereo Jack (example is an OEM parameter) Seed Adjust Seed Value Sunder Value Sunder Value Sunder Value Sunder Value Seed Seed First Stereo Sunder Value Sunder Value

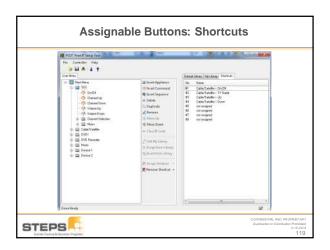

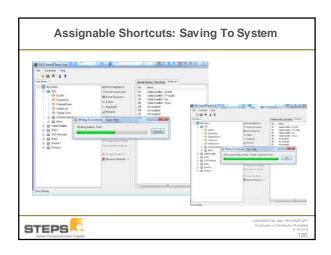





### **MWD - Pulse 6 Specifications**

- HCPCS Code
  Pulse 6 BC: K0848
  Pulse 6 CC: K0849
  Pulse 6 CC: K0856
  Drive Wheel Position
  MWD
  Speed
  6.0 mph
  Drive Wheel Tires
  13"
  Battery Type
  22 NF
  Turning Radius
  22"
  Product Weight
  235 260 lbs.

- · User Weight Capacity
  - 300 lbs.
- Overall Base Width
- 24.5" Overall Length
- 33.9"
- Seat Width
  - 12"-22"
- Seat Depth
   12" 22"
- Seat to Floor Height
  BC: 16.5",17.5", 18.5"
  CC: 17", 18", 19"
  SC: 17.5", 18", 19.5"
- STEPS

### Pulse 6: Footprint Performance Compact Footprint 33.9" STEPS

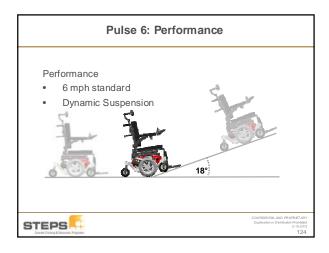

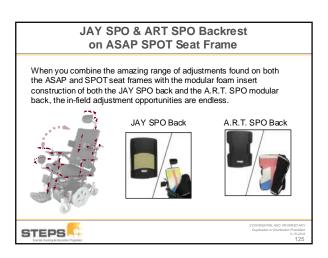



### **Pulse: Hanger Options**

- Fixed Centermount
- 90 deg
- 75 deg taper
- 70 deg Swing-away
- 65 deg Swing-away
- Power Center mount single and Ped.





CONFIDENTIAL AND PROPRIET ARY Duplication or Distribution Prohibition 3-15-2015 127

### **ASAP II Seating Adjustments**

 Check all bolts and if you run into any of the bolts not matching see bottom of Seat Frame depicted below 2 adjustment area of inside where frame inserts and outside where frame attaches.







ONFIDENTIAL AND PROPRIET AS Duplication or Distribution Prohibits 3-15-20 128

### **ASAP II Seating Adjustments**

 Check to make sure all parts are tightened and check width 16" wide complete and measures up to expectations



STEPS

CONFIDENTIAL AND PROPRIETS
Duplication or Distribution Prohib
3-15-2
12

### Cantilever Arm Adjustment

- 1. The Cantilevered Armrest angle adjust range is from 0° to 40° in 5° increments. To adjust the ammest angle up or down, losen the photo bot (A) and urscrew the 5mm both (B) using a Hex key. Adjust as necessary, using ether the inside or ouiside both fole (C) and then re-tighten bolts to specifications.
- Height adjustment for the aimrest also utilizes the 5mm Hex Key. Loosen the bolt as shown in pic 2 and move the armrest up or down.

Once desired height is acheived, re-tighten the bolts on both sides to specifications

3. For User-in-chair angle adjustments, Loosen the set screw with a 5mm Hex Key as shown in Pic 3 adjust the armrest angle and then tighten the set screw at desired angle. Pic 4.





CONFIDENTIAL AND PROPRIET AR Duplication or Distribution Prohibite 3-15-201

### **Pulse: Tilt Module**


- Tilt Module mounts between the ASAP seat and power base
- The module is replaced as an entire unit
- Tilt modules can be added after initial purchase if necessary







ONFIDENTIAL AND PROPRIETAR Duplication or Distribution Prohibits 3-15-20 13'

### 

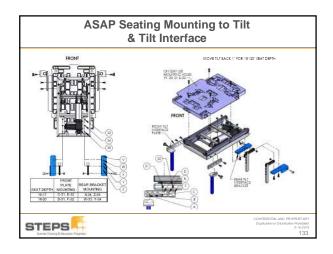

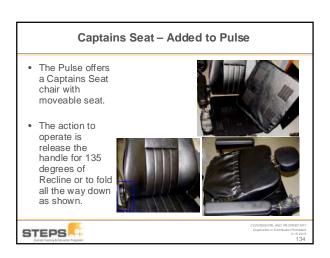

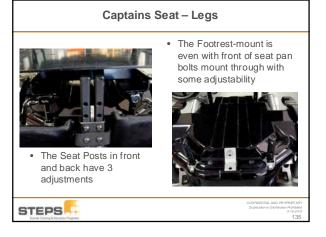

### Captains Seat - Added to Pulse

 The adjustment of bottom seat pan can be done through the top




- Or can be adjusted by moving the brackets on the side
  - It has 3" of depth adjustment in 1" increments








CONFIDENTIAL AND PROPRIET AR Duplication or Distribution Prohibits 3-15-20 4 38

### For chairs with tilt install the adaptor harness. 1. The 3-way connector plugs into the controller OBC port. 2. The 6-way connector accepts the actuator harness from the tilt. 3. The 4-way connector mates either to the TM40 flying lead or the drive through harness. Actuator Adaptor Connector VR2-90 Controller Connector View See the photos in upcoming slides See the photos in upcoming slides

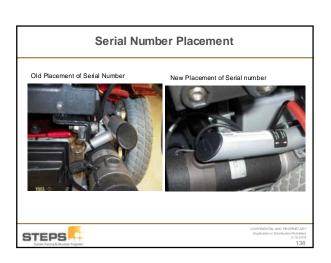

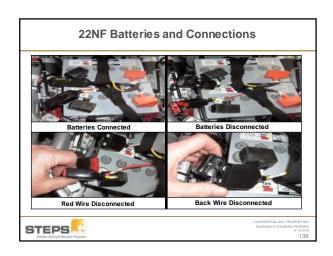

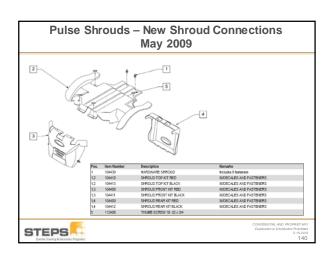

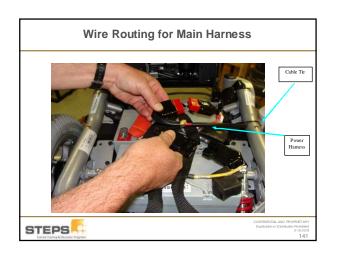

### Wire Routing for Main Harness

- Inline fuse requirement: in order to satisfy a strict reading of ANSI/RESNA section 14 and properly protect the power wiring against over-current faults, a protection device must be placed as close to the battery terminal as possible. The circuit breakers we use do not lend themselves to this, so an inline fuse near the battery post is the preferred solution. These fuses are required whether the breaker is retained or not.

  Advanced electronic controls: The motor controllers that have become industries and the leaf of the part of the leaf of the property and the leaf of the leaf of the property and the leaf of the property and the leaf of the property and the leaf of t
- Advanced electronic controls: The motor controllers that have become industrystandard over the last 20 years now handle the motor overload protection function previously handled by the circuit breaker. The onboard microprocessor tracks the motor currents over time and reduces them according to parameters determined by motor and chair testing. In fact, because of PVM technology efficiency, the circuit breaker no longer sees enough current in a stall condition to properly protect the motors. This must be done by monitoring not the battery current but the actual motor currents as the controller does.



CONFIDENTIAL AND PROPRIET ARY Duplication or Distribution Profibilities 3-15-2013

### **To Remove Motor Assembly**





Remove holt from front arm fire

Bolt with spacer & outside washe





Bolt with spacer & inside washer Don't reverse linkage small offset


STEPS

Duplication or Distribution Prohibits
3-15-20
4.41

### Pulse Motor Assembly | Table | Table

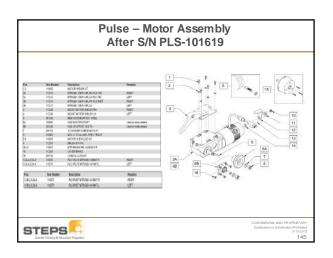

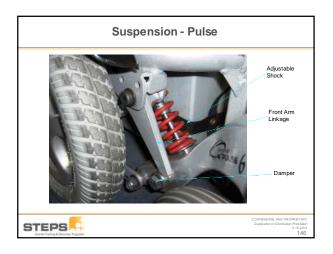

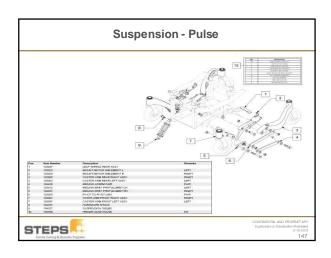







### Changes to S646 and P222SE

- Standardize Motors and Electronics across the entire power portfolio.
- Simplify overall offering.
- Increase overall reliability.
- Update electronics to a more robust solution.
- Add an LED package so Europe can access these same products. The LED package is also becoming a popular group 4 option in the US.





CONFIDENTIAL AND PROPRIET ARY Duplication or Distribution Prohibited 3-15-2015

### Overview of the Upgrade

Highlighted below are the primary changes to both the Quickie S-646 SE, Quickie S-222 SE, and in some cases the Quickie S-636

	Quickie S-636 (with standard motor package)	Quickie S-636 (with torque motor package)	Quickie \$-636 (with community motor package)	Quickie S-646 SE	Quickie P-222 SE
HCPCS Code	K0856, K0861	K0856, K0861	K0856, K0861	K0868, K0877, K0884	K0868, K0877
Speed	6.5 mph	5 mph (25%more brque)	8.5 mph	8.5 mph	8.5 mph
ser Weight Limit	300 bs	400 lbs	300 bs	250 lbs	350 lbs
Upcharge	0	\$500	\$1,000	\$0	\$0
Motors	Linix Motors	Linix Motors	Linix Motors	Linix Motors	Linix Motors
Electronics	VR2, R-net	VR2, R-net	VR2, R-net	VR2_R-net	VR2, R-net
Power Tilt.	TRAX seating with integrated Tit	TRAX seating with integrated Tilt	TRAX seating with integrated Tit	TRAX seating with integrated Tilt	Rehab Tilt
Power Recline	Standard Rehab Recline Seat	Standard Rehab Reclin e Seat	Standard Rehab Reclin e Seat	Standard Rehab Reclin e Seat	Not Available
D Light Package	Not Available	Not Available	Not Available	Available	Available
All available Quickiecolors	Available	Available	Available	Available	Available



ONFIDENTIAL AND PROPRIET AR Duplication or Distribution Prohibite 3-15-201

### **Benefits to This Change**

New Features	Benefits:	
Common Electronics Platform	Helps Reduce Inventory -Elihinate nee dib inventory Otrori de electrorics - Reduce overall in sen pry cost - Incea se familiarity	Plug and Play Capability - Carlsup R-net plastacions injust est cass who use set or reprogram chair - Omnisupp of sup 102 different input desices
Common Motor Mix	Familiarky - 7 echs lamiliar vith the same motor across multiple if net	Opens up HD option in Group 3  - The Heavy Laty option is much more pap-sier in Group 3  - The James of the American Companies in Group 3 anegory, there is an increase and potential for more business.  - Also opens a high performance option for the Outcide
Reduced Complexity	Hamesses are Modular, i.e., a Specific Harness for Each Piece of the System - Immery cases its immens the need for midgle harness is right, which have been seen to action right harnesses when moving a searing platf	Replaced Circuit Breakers with Infline Fuse poster against over hearing, but do not prosect against over hearing, but do not prosect against intigs. We a dead short C controllers not procedure gainst controllers of controllers not procedure gainst controllers of controllers not procedure gainst controllers of controllers not procedure and controllers of controller
Reliability and Performance	Proven Base Design - Over 10 years of proven performance	Proven Electronics  - No Crives VP2 and R-visity sens are so eliable to the Crives VP2 and R-visity sens are so eliable to the Crives VP2 and R-visity sens are so eliable to the Crives VP2 and VP2 a



IFIDENTIAL AND PROPRIET.

Splication or Distribution Prohib


3-15-2


### **Common Electronics Platform** Plug and Play Capability - Can swap R-net joysticks/Omni Input devices Helps Reduce Inventory - Eliminate need to inventory Otronics ${\tt Common}$ Electronics P la tfo r m More Joystick Options - More profiles - Charging through the joystick - Color display and real-time clock - IR and Bluetooth availability using the Improved Driving Algorithms - More intuitive handling - More fluid transitions STEPS Same Motor Type Familiarity -Techs familiarith the same motor across multiple lines Opens up HD option in Group 3 - The Heavy Duty option is much more popular in - rine Heavy Duty option is much more popular the Group 3 coding category. - Also opens a high performance option for the Quickle. Common Motor Mix STEPS **Reducing Complexity** Harnesses are Modular, i.e., a Specific Harness for Each Piece of the System - In many cases it diminates the ned from tliph is manse (angla, which simplifies the driving process - Reduces the need to acts mit gla harnesses when moving a seating platf Replaced Circuit Breakers with Inline Fuse - Circuit breakers practs against overheading, but on reprotect against things like a dead short Controllers nowprotect against overheading. Circuit breakers degrade some every time they are triped. Reduced Complexity STEPS

# Reliability and Proven Base Design - Over 10 years of proven performance Proven Electronics - RD Dives VE2 and R-rat systems are so reliable fait they are used by a number of prominent power wheelchair manufacturers throughout the industry Proven Motors - Less than a 15% what mate lowest in company intolly - Same grown motors assued on the OMT series, - SAD System - Same grown motors assued on the OMT series, - SAD System - Allower from Same and a Outliefe berry - After reliable break-to base design COUNTEDITIES AND MOTHER ART Displaced on the Date of Date in peach to the at security and the series of Date in the Counter of the series of Date of D



### A Power Portfolio Example: Quickie • Rear Wheel Drive (RWD) - S-636/S-646 - P222 SE • Mid Wheel Drive (MWD) - Pulse 6 - QM 7 Series

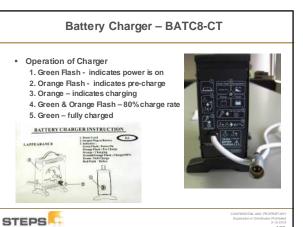


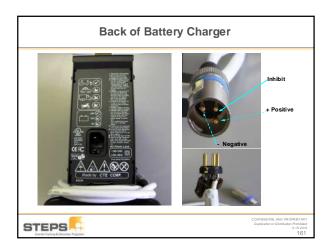







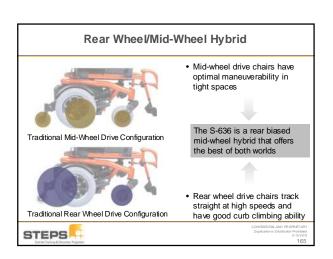

### RWD - Quickie S-636 Comparison


- Same Rear Wheel Drive Performance with Tilt OR Tilt and Recline capability as the S646
- Same Shock Suspension as the S646
- 6.5 mph ONLY w/ 0-300 lbs weight capacity
- Group 24 Batteries only
- Otherwise the same product as the S646.



STEPS

CONFIDENTIAL AND PROPRIET AF Duplication or Distribution Prohibits 3-15-20

### Lever Arm Release Old Motors Gear Release Freewheel Changing: Difficult to engage, often requires rocking to release rocking to release grease.

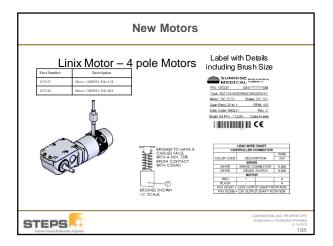







# Maintain Traction Over Transitions The dynamic anti-tip suspension can move upward to keep the rear wheels in contact with the ground during transitions. \*\*CONFERTIL AND PROPER ANY Deplication of Database Page 1, 15 200





# QM-710 - Full Featured Complex Rehab Power Chair Modular Components Available in 3 Models Quickie QM-710 (BC) – No Power Option Base with ASAP II Rehab Seat Quickie QM-710 (SC) – Single Power Option Tilt (SPOT) Base with ASAP II Rehab Seat Quickie QM-710 (MPC) – Multiple Power Option Base (SPOT Tilt and Recline) with Recline Rehab Seat

### QM-710: Motors

- Industry Proven Components Linix Motors
- Less than 1% return rate, lowest return rate in company history
- Same proven motor as used on the S-636, Pulse 6, and Rhythm
- Use unique dual end-of-line inspection system: both at vendor and at Quickie factory





STEPS

CONFIDENTIAL AND PROPRIET A Duplication or Distribution Prohibi 3-15-21

### • Shroud Available in 5 Colors \*\*CONFERRATION AND PROPRETARY Digital along the National Property and the National Propert



### QM-710: Spider Track Suspension

- Gas springs provide a dampen to reduce the jolt of an impact
- Example of a spring suspension
- Example of a dampened suspension

Pure Spring Suspension Very Bouncy



Classic Cadillac Dampened Suspension Very Smooth





CONFIDENTIAL AND PROPRIET AS Duplication or Distribution Prohibit 3-15-26 17

### QM-710: Spider Track Suspension

• When you slowly lower in a pneumatic office chair, it is the gas release that dampens the fall and/or impact.

Pneumatic Office Chair – classic example of a dampening system





STEPS

ONFIDENTIAL AND PROPRIE Duplication or Distribution Prof 3-15

### Triggered Gas Spring Dampening

- PREVENTS collapse of the front casters
- PREVENTS shock associated with non-suspension fixed front casters





CONFIDENTIAL AND PROPRIET AR Duplication or Distribution Prohibite 3-15-201

### **Triggered Gas Spring Dampening**

- Gas dampening system, prevents collapse but dampens the impact
- User maintains seated position during transitions
  - Dampened contact reduces tonal responses
  - Helps prevent user displacement
  - Helps ensure a stable seating posture



STEPS

CONFIDENTIAL AND PROPRIET A Duplication or Distribution Prohib 3-15-2 18

### QM-710: Traction & Articulation

- Climbs 3" curbs with a 22° transition angle
  - Design maximizes front caster articulation range








STEPS

CONFIDENTIAL AND PROPRIE Duplication or Distribution Prol 3-15

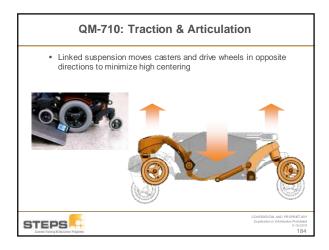

### QM-710: Designed With Minimal Parts Complexity

- 33% reduction in parts compared to legacy Quickie high-end power chairs
- Fewer parts reduces
  - Fewer parts to potentially fall out of tolerance, which reduces potential breakdowns
  - · No cable adjustments
  - Fewer parts to service creates a reduced service turnaround
  - · Easier for technicians to understand





ONFIDENTIAL AND PROPRIET A Duplication or Distribution Prohib 3-15-2 18

### SpiderTrac Suspension: Gas Strut

- Suspension mechanism accomplishes two tasks with one mechanism
- Suspension <u>prevents</u> forward pitching and high centering with a single mechanism



The rear casters rotate downward as the chair goes over a curb, this triggers a constriction of gas flow in the gas spring, which dampens front caster impact and prevents caster collapse.

The linkage between the front caster and rear caster prevent high-centering. The same gas spring is used to absorb shock for both casters.





CONFIDENTIAL AND PROPRIET A Duplication or Distribution Prohibi 3-15-21

### **Removal of Pneumatic Gas Piston**

- Gas Spring Removal

  1. Remove the wheel assembly by using a 5/8\* deep socket wemen to remove the 3 bolts that hold the wheel and inn in place (figure1). Before removing the Gas Spring raise the certer of the chair with a
- For Gas Spring Removal, Remove the Lower holding bolt (A), using a 5mm Hex key (figure 2), Note: The shocks used in this power base are pre-adjusted at the factory.
- Remove the top holding Bolt with a 6mm hex in Pic. 3.







### **Removal of Pneumatic Gas Piston**


- 4. Once the top holding bolt is removed, loosen the set screw (A) as shown by using a 4mm Hex wrench (figure 4)
- 5. Once set screw is loosened, the gas spring shaft can be unscrewed and removed from the securing point as shown in (figure 5)
- 6. To re-install the Gas Spring, Reverse the or to enistant use Gas spining, reverse the gas prings threaded shaft to the correct distance of 12-13mm for required installation clearance (C) as shown in (figure 6). Then tighten the Jam Nut (D) to secure the shaft at the correct distance, and replace the gas spring insert.










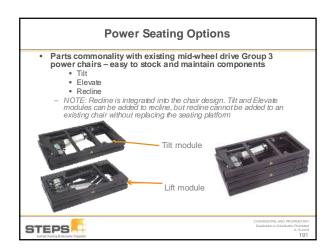

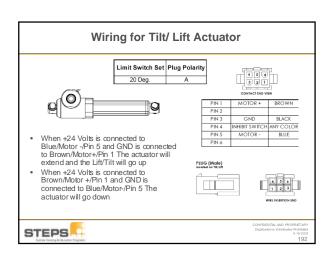



### 3 Seating Module Choices

### Non-expandable Electronic Option for Challenging Funding

- Non-expandable 3-Axis control available in three options
  - 2-Axis tilt and lift (direct activation)
  - 3-Axis tilt and power legrests (toggle activation)
  - 3-Axis lift and power legrests (toggle activation)

2-Axis box 4.5 & 7.0 Amp box used for Lift and Tilt 3-Axis 7.0 Amp box us ed for Tilt and Pwr Legrests 3-Axis 4.5 Amp box us ed for Lift and Pwr Legres ts









CONFIDENTIAL AND PROPRIET A Duplication or Distribution Prohibit 3-15-20 19

### **Diagnostics Overview**

- 1. Saving files to send via E-mail etc to make product work again.
- 2. Saving files and Diagnostic Log and opening on P.C. look then load back to chair.
- 3. Understanding Inhibits what is causing that flashing mode or Turtle on Screen for LCD models.
   4. Deal from Models from models.
- 4. Real time Monitoring of Motor Current Motor Voltage Inhibit Operation in Troubleshooting – including Seating controls and battery voltage.
- 5. Joystick Monitoring out of neutral switch activation.
- 6. Reading proper ohms resistance know where and when to check resistances.
- Resources such as Tech Support On-Line center for Videos –
  Setup Guides OBP programming Road Maps etc.



ON FIDENTIAL AND PROPRIET A Suplication or Distribution Prohibi 3-15-21

### **Diagnostic Procedure**

- · Please follow this procedure:
  - Read and note the Trip Text displayed, the identified Module and the Trip Code. Switch off the control system.
- Make sure that all connectors on the listed Module and the wheelchair are mated securely.
- Check the condition of the battery.
- Find the definition of the Trip Code in this Service Guide, or go online to PGDT.com and insert code and take required action.
- Switch on the control system again and try to drive the wheelchair. If the safety circuits operate again, switch off and do not try to use the wheelchair. Contact Sunrise Technical Service.



CONFIDENTIAL AND PROPRIET ARY Duplication or Distribution Prohibited 3-15-2013

### Display & Screen Diagnostics



- Display shows the details that are needed for a Technician to troubleshoot the problems.
- Identified Module Identifies which module has the present problem
  - PM = Power Module
  - JSM = Joystick Module
  - ISM Intelligent Seating Module
  - OMNI Specialty Control Device and Display
  - ATT attendant Control



CONFIDENTIAL AND PROPRIET AS Duplication or Distribution Prohibit 3-15-26 19

### **Display & Screen Diagnostics**

### Trip Text

Text Version of the code Trip Code



- Is the Alphanumeric Code indicating the designated code number letter combination.
- In the display we have a Power Module with a Trip Text of Low Battery and a Trip Code of 2C00 – depending on type of programming device it will store either a combination of trip code and Text or Text only.
- Programmers prior only displayed Trip Codes sometimes making it hard to decide where and what to look for.



CONFIDENTIAL AND PROPRIET AR Duplication or Distribution Prohibite 3-15-201

### **VR2 Diagnostics**

### Using PP1 or Qtronix Programmer

- · Read System Log?
- The control systems have a diagnostics log facility which stores the number of occurrences of the last eight trip codes. This allows you to view the contents. The display format is as below.
  - 1: Code 2C00, #1
  - 2: Code 3C00, #3
- No more entries
- This reads line by line as.
  - Line 1 trip code 2C00 has occurred once
  - Line 2 trip code 3C00 has occurred three times

LED Diagnostics - R-Net EL Control

· Only two trip types recorded.



0

Φ

Ф

### LED Control System Status Indication: The battery gauge and Maximum Speed/Profile indicator show the status of the control System. With the help of this listing it is possible to diagnose problems that may appear drastic at first, but may have a very simple and easy to fix cause. 1. Battery Gauge is Stead y. This indicates all is well Battery Gauge Flashes Slowly. The control system is functioning correctly, but you should charge the battery as soon as possible. 3. BatteryGauge Steps up: The wheelchair batteries are being charged you will not be able to drive until chair is disconnected from charger. Switched Power On and Off again. Battery Gauge Flashes Rapidly (even when joystick released): The control system safety circuits have operated and the control system has been prevented from moving the chair. This indicase a system trip. The R-N-the has deleted a problem somewhere in the wheelchairs electrical system. Please follow below procedure: 1. Switch of the cortrol system. Selb-1 of the contral system. Sees when we show it is colour. 2. Make are that all connect or so in the wheelthat and the contral system are made correctly. 3. Check the contributed the battery. 4. Seach the contributed many and try to drive the wheelthat. If the safety circuits oper all again. 5. Speed Indicator Ripples Outwards From Center: In this instance the LEDS make a prigher motion starting with the Middle LED and then stepping outwards on both sides. The Control System has descend that a new module has been added and is re-configuring. Please allow for a few seconds for this to complete.


- 0


Speed Indicator LED's 2.8.4 Flash: Then the control system requires a reboot. For example, after a module re-configuration, the second and fourth speed indicator LED's will

STEPS

### LED Diagnostics - R-Net EL Control Note: If the Programmable Parameter, Motor Swap has been enabled, the left and right hand references in the table will be 7 Bar A joydox tour is indicated. Wasy sure that the popular is in the center position before self-control partiers. the polaring brakes have a bad connection, check the passing brake and motor or use the control system connections are secure. An excessive votage has oven applied to the control system. This is usually of connection. Check the bottopy connections. A communication tout is indicated. More sure that journey cases is securely connection age? STEPS

### Indication of Display Related Faults R-net LED Joystick 1 LED The battery needs charging or there is a bad connection to the battery. Check the connections to the battery. If connections are good, try charging the battery. 2 LED The M1 Motor – has a bad connection. Check connections and resistance of motor. 3 LED The M1 Motor – has a short circuit to a battery connection. Motor short can lead to damage of control module. 4 LED The M2 Motor – has a bad connection. Check connections and resistance of motor.

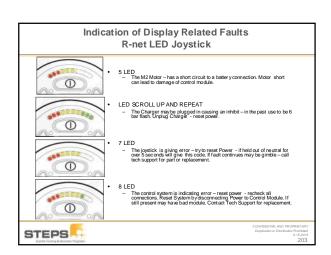

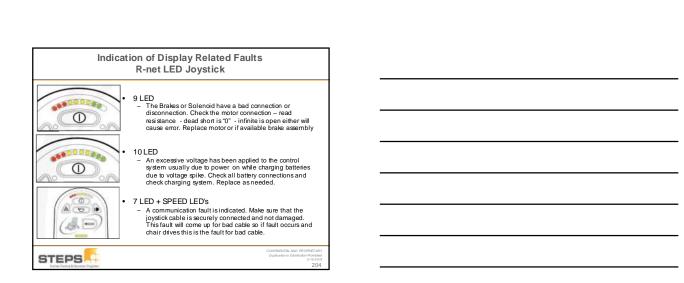

### Indication of Display Related Faults R-net LED Joystick

Actuator LEDs Flash
 An Actuator Trip is indicated. If more than one actuator is attached, check which actuator is not working correctly. Check for resistance and if any thing is causing excessive current.



- Maximum Speed / Profile Indicator Flashes
   This indicates the speed of the wheelchair is being limited for safety reasons. The exact reason will depend on if a seating control is on the chair or not. The most common cause is seat is elevated position.
   Maximum Speed / Profile indicator Ripples Up and Down

  This indicate the control output is before. Pofer to pay the page for unlocking.
- - This indicates the control system is locked. Refer to next page for unlocking procedure.



### **Diagnostic for Unlocking Unit**



If an LED Joyatek Module is use (EL) the Speed Indicator LED's will ripple from left to right
To unifor the wheelchair:

If the control system has switched off, press the Outoff button.

Deflect the joyatic freewards until the control system beeps.

Deflect the joyatic freewards until the control system beeps.

Deflect the joyatic freeward is a long beep.

Release the joyatic fire will be a long beep.

The wheelchair is now unlocked and ready for use



### **Diagnostic Text Definitions**

- The R-net has various connections that are referred to by the diagnostic text.

  Power Module (PM)

  Joystick Module (JSM)

  - Intelligent Seating/Lighting Module (ISM)
















### Thank You For Attending! Ken Kalinowski, Senior Service Technician ken.kalinowski @sunmed.com


STEPS

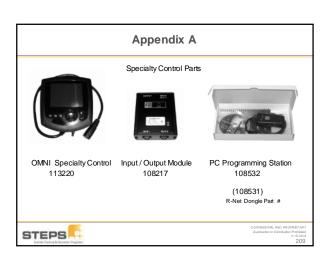




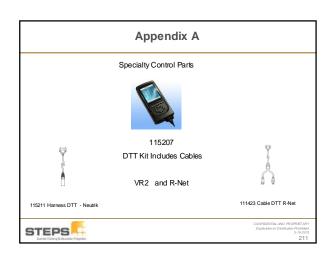


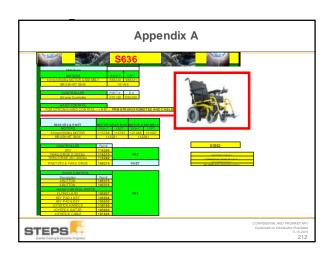

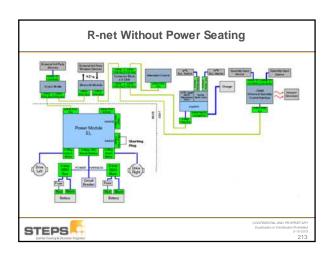

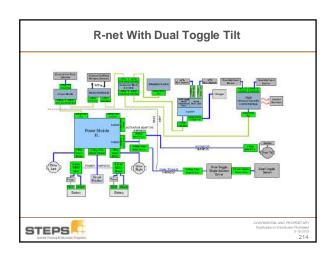



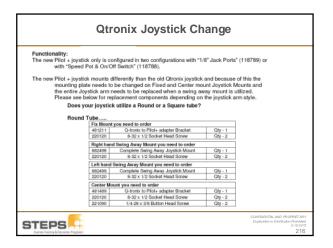

	Square Tube									
		vou need to orde	ır							
	481211				Oty - 1					
	220120		Socket Head Screw Oty - 2							
	Dight han	d Swing Away Mo	e and work man	ad to order						
	118824	Complete Swir			Oty - 1					
	220120		Socket Head		Oty - 2					
	118825	Left hand Swing Away Mount you need to order  118825 Complete Swing Away Joyntick Mount City - 1			Oty - 1					
	220120				Oty - 2					
	Center Mount you need to orde 481489 O-tronk to Pilot			e Dranket	Oty - 1					
	220120		Socket Head Screw Oty - 2							
	221090		Button Hea		Oty - 2					
	-									
			Plots PN				Feature Comparison			
Qtronix PN	Des criptic	Description		D	escription	Pot	Toggle	Jacks	Right	Left
									Hand	Hand
139515	REMOTE BOX W/O JACK R	OTE BOX W/O JACK R S525 11		JOYSTICK P+ PC	T/TOG NO JACKS	- 4	4	NA	- 4	NA
139516	REMOTE BOX W/JACKR S	REMOTE BOX WUACKR \$525		JOYSTICK P+ JA	CKSNO POT/TOG	X		- V	- 4	NA
139517	REMOTE BOX W/O JACK L SS25		118788	JOYSTICK P+ POT/TOG NO JACKS		4	4	NA	NA	х
139518	REMOTE BOX W/JACKL S525		118789	JOYSTICK P+ JACKSNO POT/TOG		×		V	NA.	х
							_			
						4	Feature R			
						X	Feature I	ost		
						NA.	F		on Otronia	

